


ADFS852

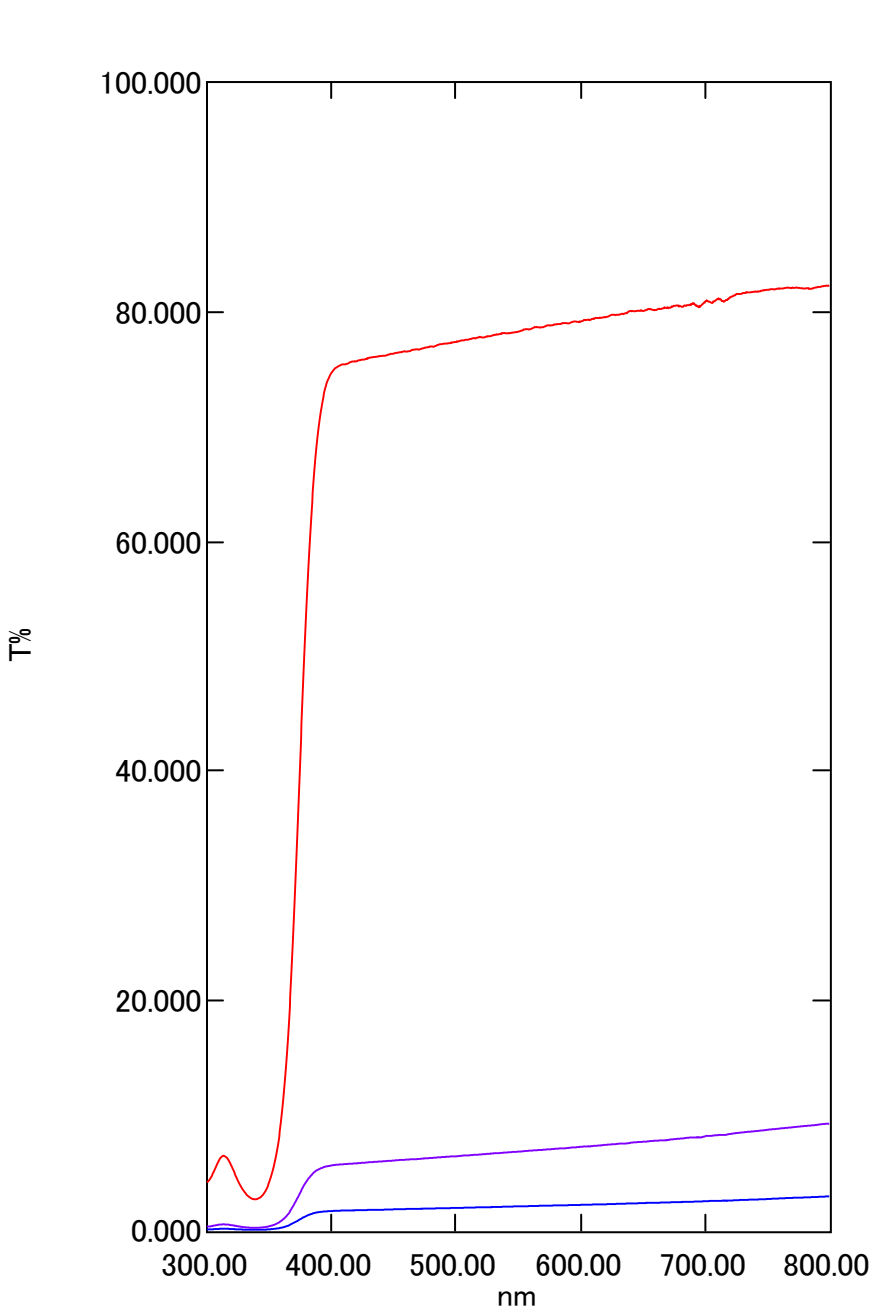

— ADFS852-0 - 0mm
— ADFS852-100 - 100mm
— ADFS852-50 - 50mm

波長 (nm)	0mm ...	50mm ...	100mm ...
300.00	4.209	0.347	0.107
310.00	6.200	0.527	0.157
320.00	5.537	0.471	0.139
330.00	3.376	0.298	0.090
340.00	2.747	0.247	0.074
350.00	4.211	0.383	0.114
360.00	10.014	0.874	0.263
370.00	27.567	2.285	0.689
380.00	54.502	4.300	1.291
390.00	70.215	5.356	1.612
400.00	74.701	5.674	1.710
410.00	75.425	5.788	1.749
420.00	75.700	5.858	1.775
430.00	75.974	5.955	1.801
440.00	76.139	6.026	1.824
450.00	76.352	6.105	1.854
460.00	76.528	6.176	1.882
470.00	76.712	6.238	1.902
480.00	76.977	6.332	1.927
490.00	77.212	6.413	1.952
500.00	77.369	6.470	1.975
510.00	77.566	6.553	2.002
520.00	77.790	6.635	2.030
530.00	77.965	6.725	2.052
540.00	78.134	6.805	2.082
550.00	78.259	6.878	2.111
560.00	78.498	6.967	2.139
570.00	78.664	7.039	2.167
580.00	78.904	7.123	2.194
590.00	78.977	7.196	2.216
600.00	79.128	7.294	2.242
610.00	79.377	7.369	2.272
620.00	79.521	7.466	2.296
630.00	79.741	7.562	2.337
640.00	80.076	7.640	2.376
650.00	80.068	7.738	2.388
660.00	80.147	7.823	2.430
670.00	80.367	7.894	2.453
680.00	80.517	8.008	2.488
690.00	80.698	8.114	2.516
700.00	80.810	8.199	2.559
710.00	81.104	8.316	2.595
720.00	81.238	8.419	2.622
730.00	81.583	8.547	2.669
740.00	81.748	8.655	2.707
750.00	81.910	8.776	2.751
760.00	82.001	8.888	2.804
770.00	82.085	8.995	2.852
780.00	82.022	9.108	2.887
790.00	82.122	9.207	2.936
800.00	82.298	9.317	2.986

[測定プロパティ]
 波長範囲(nm): 300.00 ~ 800.00
 スキャン速度: 高速
 サンプリングピッチ: 0.5
 オート サンプリングピッチ: 有効
 測定モード: オート

[装置プロパティ]
 装置の種類: UV-3600Plus シリーズ
 測光値: 透過率
 スリット幅: (20)
 時定数: 0.1 秒
 光源: 自動
 光源切替波長: 290.00 nm

.....