


ND30AH

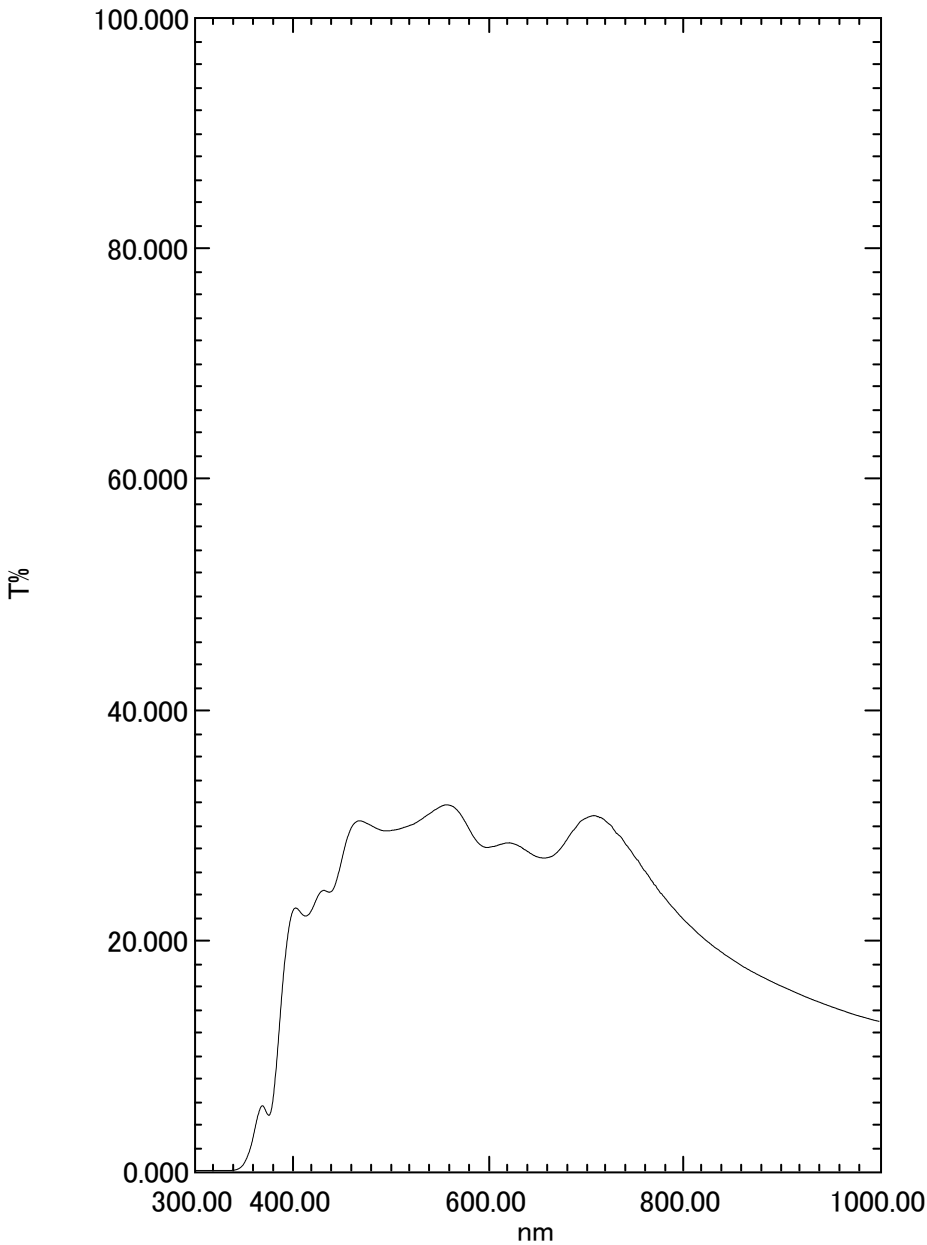

[測定プロパティ]

波長範囲(nm): 300.00 ~ 1000.00  
 スキャン スピード: 中速  
 サンプルングピッチ: 0.5  
 オート サンプルングピッチ: 有効  
 測定モード: オート

[装置プロパティ]

装置の種類: UV-3600Plus シリーズ  
 測光値: 透過率  
 スリット幅: (20)  
 時定数: 0.1 秒  
 光源: 自動  
 光源切替波長: 282.00 nm

.....

波長 (nm)	RawData ...
300.00	-0.000
310.00	0.000
320.00	0.001
330.00	0.001
340.00	0.031
350.00	0.601
360.00	3.126
370.00	5.557
380.00	6.339
390.00	16.510
400.00	22.499
410.00	22.206
420.00	22.719
430.00	24.264
440.00	24.264
450.00	26.953
460.00	29.765
470.00	30.329
480.00	29.937
490.00	29.539
500.00	29.490
510.00	29.663
520.00	29.949
530.00	30.387
540.00	30.970
550.00	31.526
560.00	31.697
570.00	31.075
580.00	29.676
590.00	28.391
600.00	28.028
610.00	28.224
620.00	28.419
630.00	28.212
640.00	27.676
650.00	27.213
660.00	27.144
670.00	27.593
680.00	28.670
690.00	29.829
700.00	30.581
710.00	30.768
720.00	30.311
730.00	29.379
740.00	28.366
750.00	27.206
760.00	26.009
770.00	24.788
780.00	23.713
790.00	22.690
800.00	21.750
810.00	20.938
820.00	20.158
830.00	19.467
840.00	18.850
850.00	18.266
860.00	17.725
870.00	17.243
880.00	16.804
890.00	16.377
900.00	15.982
910.00	15.606
920.00	15.234
930.00	14.878
940.00	14.555
950.00	14.237
960.00	13.938
970.00	13.649
980.00	13.384
990.00	13.144
1000.00	12.946