


ND50AH

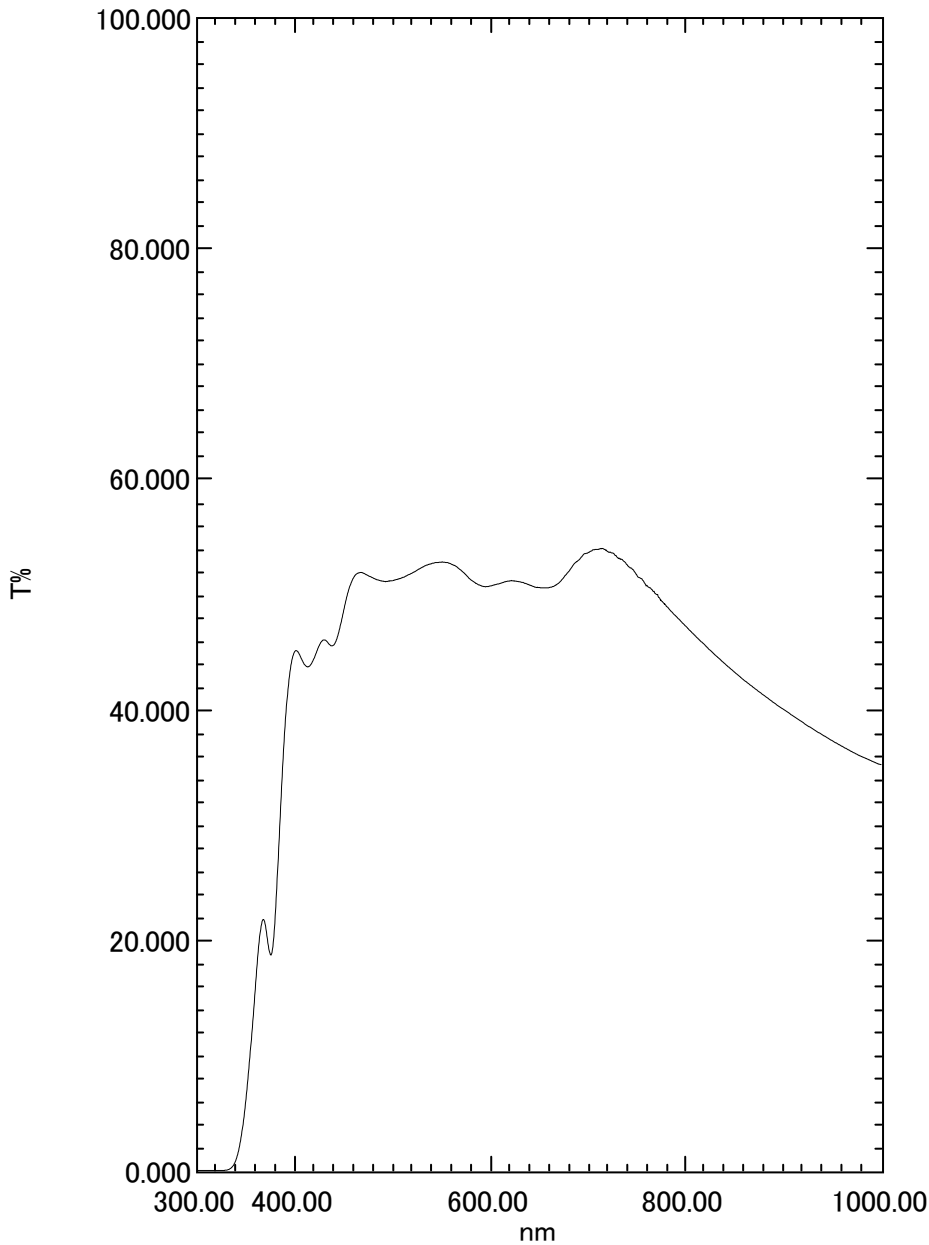

[測定プロパティ]

波長範囲(nm): 300.00 ~ 1000.00  
 スキャン スピード: 中速  
 サンプルングピッチ: 0.5  
 オート サンプルングピッチ: 有効  
 測定モード: オート

[装置プロパティ]

装置の種類: UV-3600Plus シリーズ  
 測光値: 透過率  
 スリット幅: (20)  
 時定数: 0.1 秒  
 光源: 自動  
 光源切替波長: 282.00 nm

.....

波長 (nm)	RawData ...
300.00	0.000
310.00	0.001
320.00	0.001
330.00	0.033
340.00	0.973
350.00	6.189
360.00	16.296
370.00	21.219
380.00	22.113
390.00	38.565
400.00	45.032
410.00	43.899
420.00	44.508
430.00	46.034
440.00	45.616
450.00	48.621
460.00	51.425
470.00	51.833
480.00	51.418
490.00	51.136
500.00	51.178
510.00	51.446
520.00	51.831
530.00	52.314
540.00	52.666
550.00	52.784
560.00	52.611
570.00	52.071
580.00	51.253
590.00	50.752
600.00	50.716
610.00	50.932
620.00	51.157
630.00	51.082
640.00	50.803
650.00	50.591
660.00	50.588
670.00	50.978
680.00	51.965
690.00	52.944
700.00	53.618
710.00	53.876
720.00	53.695
730.00	53.198
740.00	52.545
750.00	51.645
760.00	50.791
770.00	49.962
780.00	48.937
790.00	48.075
800.00	47.213
810.00	46.350
820.00	45.539
830.00	44.725
840.00	43.959
850.00	43.227
860.00	42.502
870.00	41.856
880.00	41.214
890.00	40.585
900.00	39.996
910.00	39.420
920.00	38.850
930.00	38.297
940.00	37.784
950.00	37.270
960.00	36.793
970.00	36.328
980.00	35.904
990.00	35.529
1000.00	35.211