


ND60AH

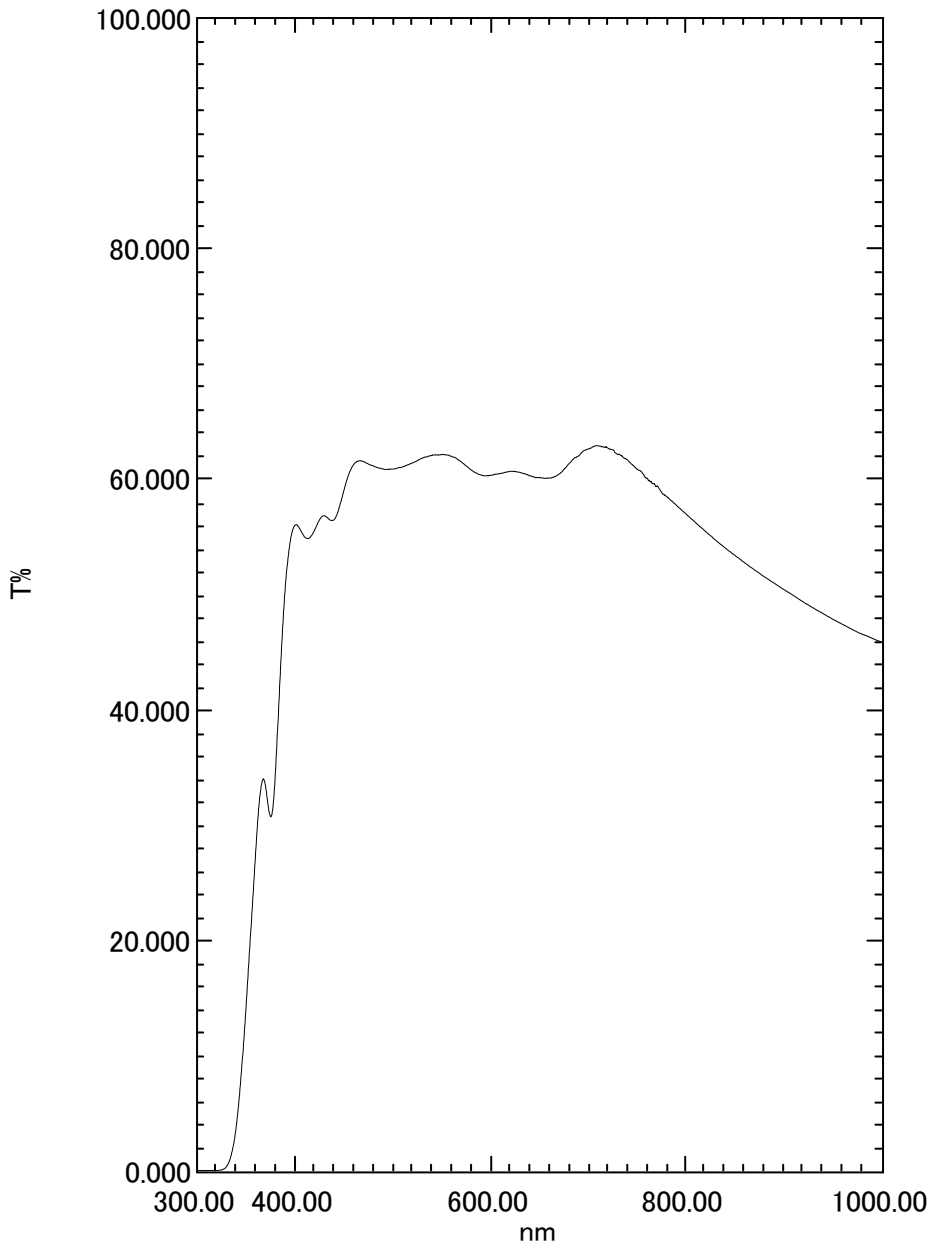

[測定プロパティ]

波長範囲(nm): 300.00 ~ 1000.00
 スキャン スピード: 中速
 サンプルングピッチ: 0.5
 オート サンプルングピッチ: 有効
 測定モード: オート

[装置プロパティ]

装置の種類: UV-3600Plus シリーズ
 測光値: 透過率
 スリット幅: (20)
 時定数: 0.1 秒
 光源: 自動
 光源切替波長: 282.00 nm

.....

波長 (nm)	RawData ...
300.00	0.001
310.00	0.001
320.00	0.005
330.00	0.328
340.00	3.790
350.00	14.093
360.00	27.693
370.00	33.374
380.00	34.286
390.00	50.304
400.00	55.968
410.00	54.984
420.00	55.542
430.00	56.804
440.00	56.464
450.00	58.957
460.00	61.238
470.00	61.521
480.00	61.142
490.00	60.879
500.00	60.855
510.00	61.052
520.00	61.388
530.00	61.792
540.00	62.054
550.00	62.123
560.00	61.963
570.00	61.505
580.00	60.798
590.00	60.360
600.00	60.311
610.00	60.473
620.00	60.649
630.00	60.577
640.00	60.283
650.00	60.119
660.00	60.076
670.00	60.406
680.00	61.295
690.00	61.998
700.00	62.599
710.00	62.877
720.00	62.688
730.00	62.130
740.00	61.679
750.00	60.841
760.00	60.075
770.00	59.386
780.00	58.472
790.00	57.746
800.00	56.969
810.00	56.206
820.00	55.447
830.00	54.715
840.00	54.033
850.00	53.380
860.00	52.737
870.00	52.123
880.00	51.531
890.00	50.965
900.00	50.419
910.00	49.884
920.00	49.340
930.00	48.820
940.00	48.334
950.00	47.847
960.00	47.399
970.00	46.950
980.00	46.546
990.00	46.197
1000.00	45.881