


ND70AH

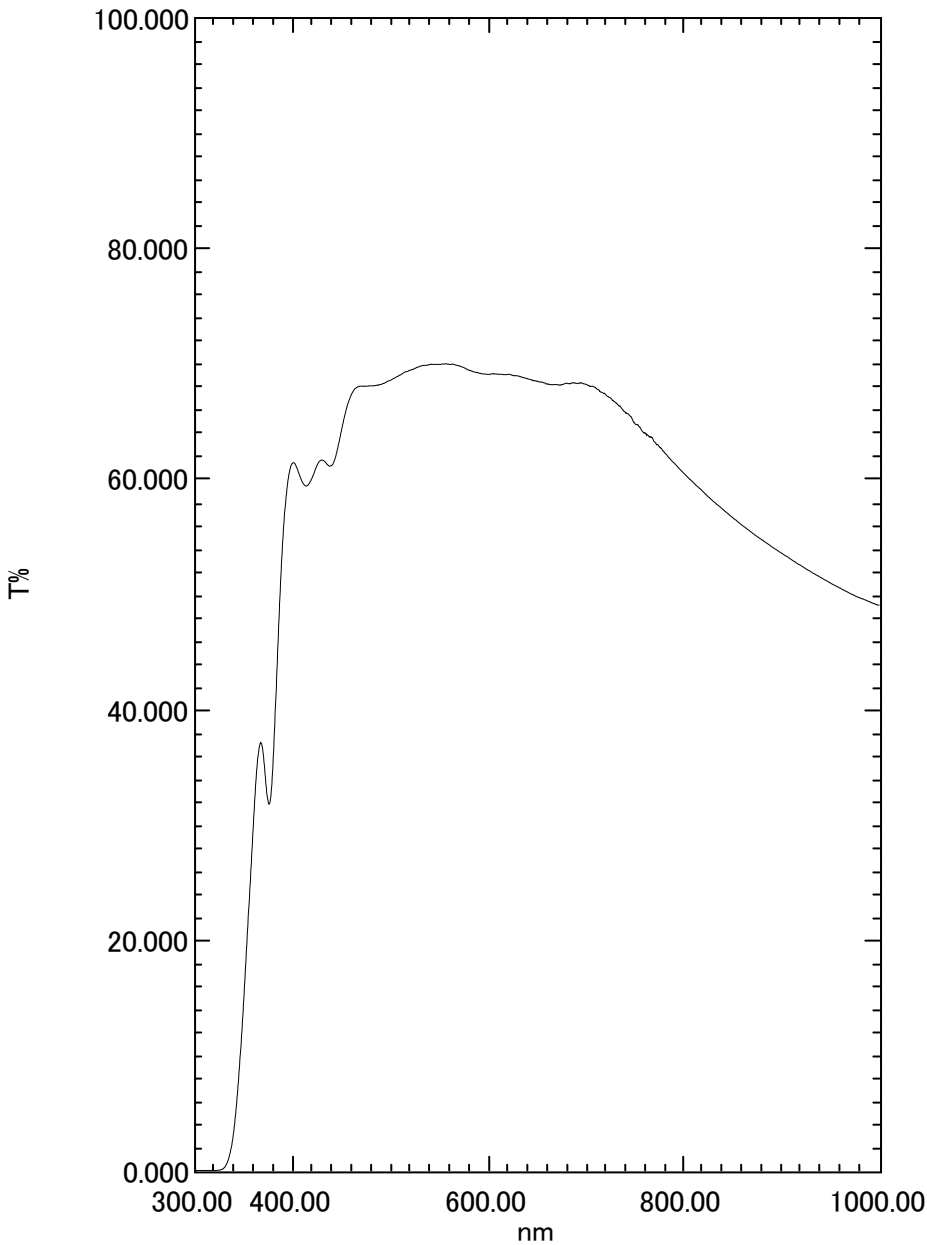

[測定プロパティ]

波長範囲(nm): 300.00 ~ 1000.00
 スキャン スピード: 中速
 サンプルングピッチ: 0.5
 オート サンプルングピッチ: 有効
 測定モード: オート

[装置プロパティ]

装置の種類: UV-3600Plus シリーズ
 測光値: 透過率
 スリット幅: (20)
 時定数: 0.1 秒
 光源: 自動
 光源切替波長: 282.00 nm

波長 (nm)	RawData ...
300.00	0.002
310.00	0.001
320.00	0.002
330.00	0.263
340.00	3.727
350.00	15.150
360.00	30.604
370.00	35.934
380.00	35.765
390.00	55.048
400.00	61.401
410.00	59.703
420.00	60.157
430.00	61.643
440.00	61.196
450.00	64.328
460.00	67.351
470.00	68.049
480.00	68.085
490.00	68.204
500.00	68.559
510.00	69.035
520.00	69.431
530.00	69.779
540.00	69.938
550.00	69.949
560.00	69.952
570.00	69.813
580.00	69.445
590.00	69.205
600.00	69.082
610.00	69.084
620.00	69.086
630.00	68.956
640.00	68.690
650.00	68.496
660.00	68.249
670.00	68.201
680.00	68.303
690.00	68.294
700.00	68.226
710.00	67.897
720.00	67.323
730.00	66.640
740.00	65.797
750.00	64.809
760.00	64.025
770.00	63.223
780.00	62.284
790.00	61.377
800.00	60.505
810.00	59.661
820.00	58.864
830.00	58.068
840.00	57.360
850.00	56.641
860.00	55.959
870.00	55.315
880.00	54.708
890.00	54.118
900.00	53.554
910.00	53.012
920.00	52.461
930.00	51.942
940.00	51.452
950.00	50.965
960.00	50.514
970.00	50.074
980.00	49.686
990.00	49.343
1000.00	49.036